



SINGLE USE



S6000

**Six Piece |
Restorative Kit**

BB56 | Ball Burnisher

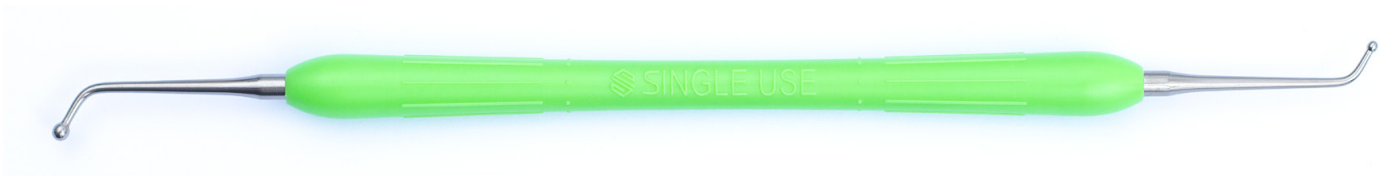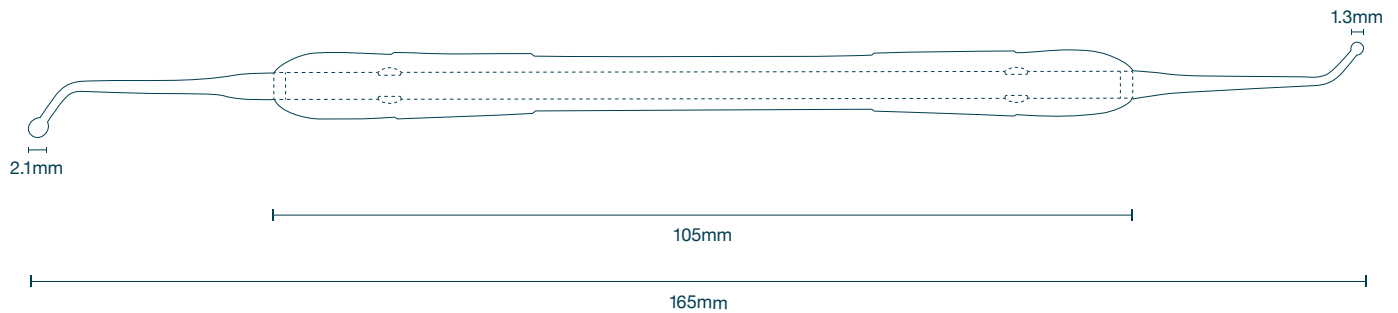

Description


The perfect blending tool. Our ball burnisher has both small and large tip ends to effectively condense and shape the composite against the cavity wall surfaces.

Features

- Double-ended instrument
- Precision manufactured stainless steel tips
- Ergonomic handle with ridges and grooves for improved grip
- Stainless steel stem for added stability
- Latex-free
- Arrives sterile and ready for use
- GS1 Data matrix
- CE Certified

Materials

- Stainless steel
- Acrylonitrile butadiene styrene (ABS)

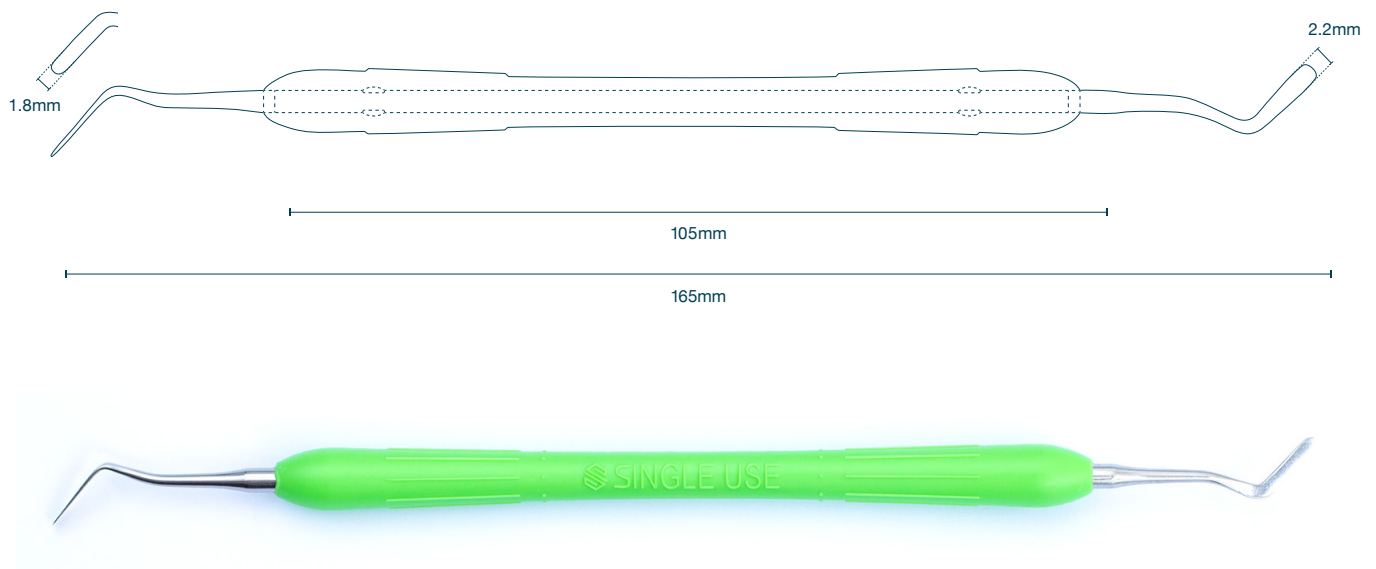

Description

Our double ended Flat Plastic features long, flat tips that are used to carry and shape composite materials.


Features

- Double-ended instrument
- Precision manufactured stainless steel tips
- Tip diameters; 2.2mm/1.8mm
- Ergonomic handle with ridges and grooves for improved grip
- Stainless steel stem for added stability
- Latex-free
- Arrives sterile and ready for use
- GS1 Data matrix
- CE Certified

Materials

- Stainless steel
- Acrylonitrile butadiene styrene (ABS)

XC67 | Small Spoon Excavator

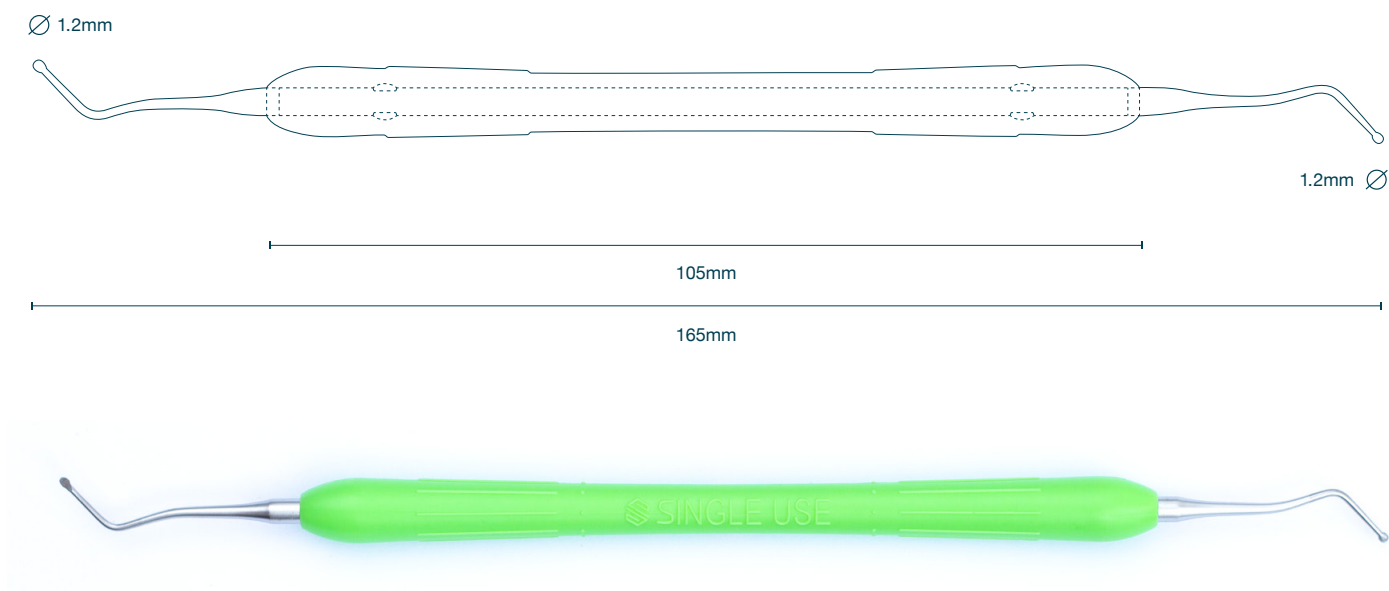

Description

One of our bestsellers, our Small Spoon Excavator is a popular choice for removing caries within deep cavities.


Features

- Double-ended instrument
- Precision manufactured stainless steel tips
- Small spoon diameter; 1.2mm
- Ergonomic handle with ridges and grooves for improved grip
- Stainless steel stem for added stability
- Latex-free
- Arrives sterile and ready for use
- GS1 Data matrix
- CE Certified

Materials

- Stainless steel
- Acrylonitrile butadiene styrene (ABS)

XC58 | Medium Spoon Excavator

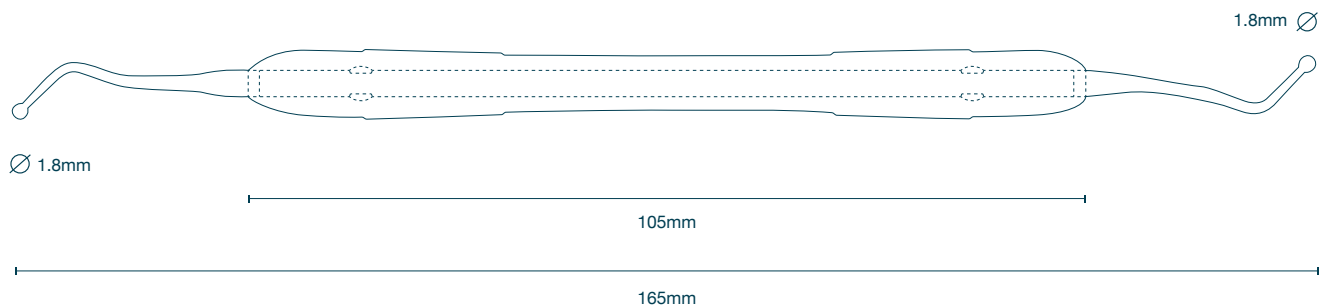

Description

Gone are the days of dull and blunt instruments. Our single-use XC58 will have sharp edges each and every time. Precision made, our medium spoon excavator adds convenience to your procedures.


Features

- Double-ended instrument
- Precision manufactured stainless steel tips
- Ergonomic handle with ridges and grooves for improved grip
- Stainless steel stem for added stability
- Latex-free
- Arrives sterile and ready for use

Materials

- Stainless steel
- Acrylonitrile butadiene styrene (ABS)

CV66 | 1/2 Hollenback

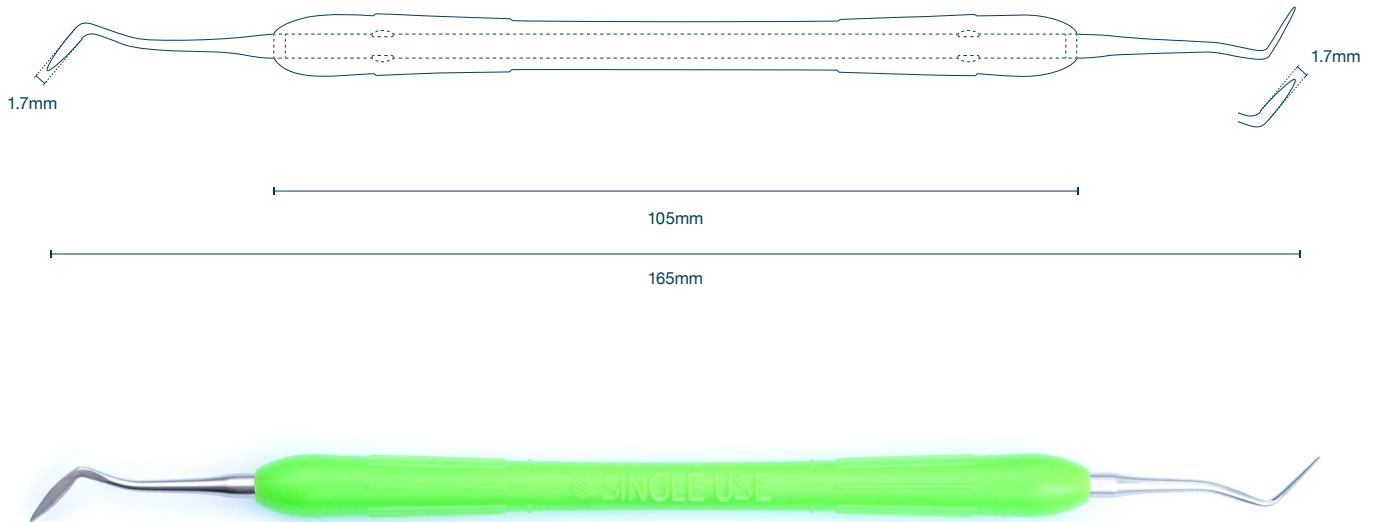

Description

Our 1/2 Hollenback is a versatile instrument used to place, carve and contour amalgam.


Features

- Double-ended instrument
- Precision manufactured stainless steel tips
- Ergonomic handle with ridges and grooves for improved grip
- Stainless steel stem for added stability
- Latex-free
- Arrives sterile and ready for use
- GS1 Data matrix
- CE Certified

Materials

- Stainless steel
- Acrylonitrile butadiene styrene (ABS)

BB49 | PFI49 Burnisher

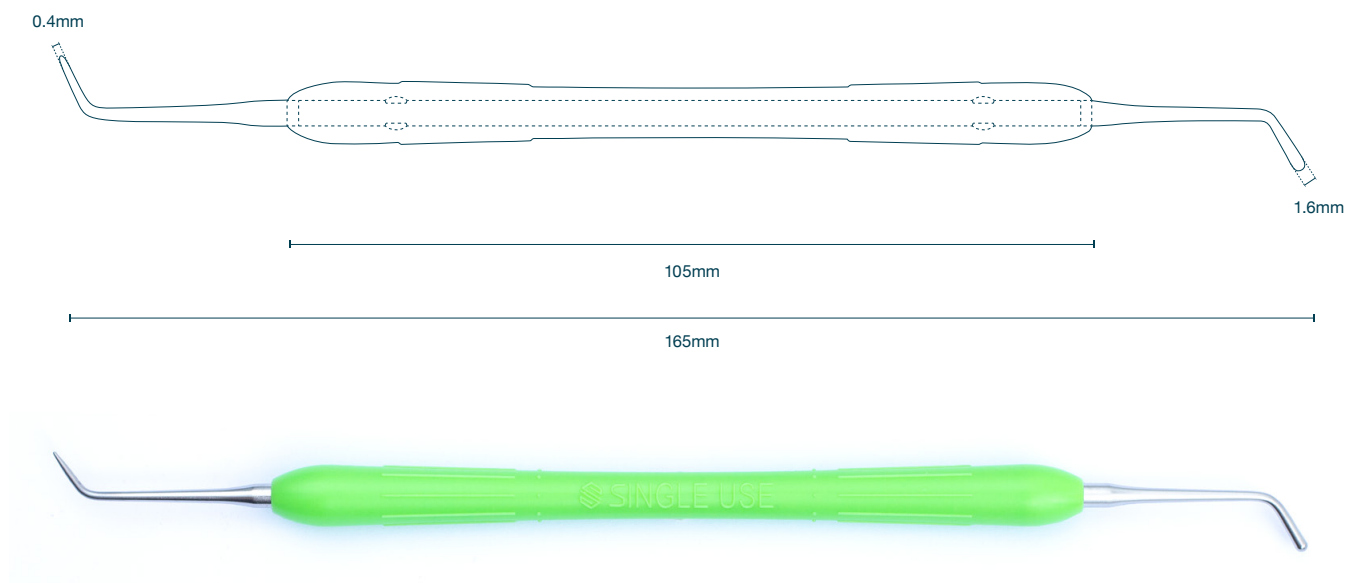

Description

The PFI49 burnisher is ideal for placing and contouring restorations. Having both a rounded and pointed tip makes this instrument useful for placing composite in small increments prior to curing.

Features


- Double-ended instrument
- Precision manufactured stainless steel tips
- Ergonomic handle with ridges and grooves for improved grip
- Stainless steel stem for added stability
- Latex-free
- Arrives sterile and ready for use
- GS1 Data matrix
- CE Certified

Materials

- Stainless steel
- Acrylonitrile butadiene styrene (ABS)

Products arrive in a protective, sterile sleeve allowing for immediate use.

Removable label with GS1 Data Matrix and product code for electronic data recording and patient logbooks.

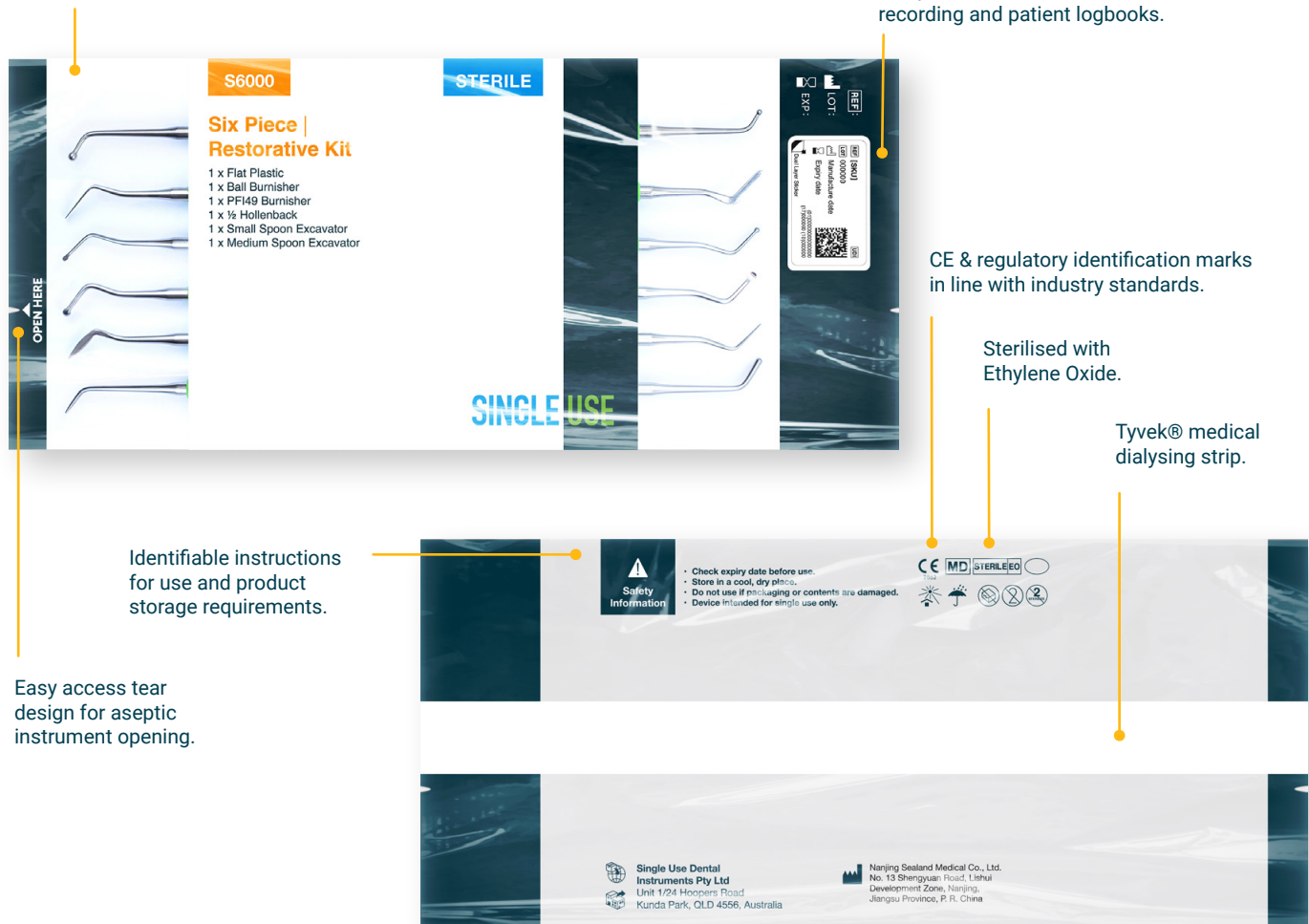

Identifiable instructions for use and product storage requirements.

Sterilised with Ethylene Oxide.

Tyvek® medical dialysing strip.

Easy access tear design for aseptic instrument opening.

Our Products

Our unique, single-use instruments were designed with simplicity in mind. Featuring a broad variety of diagnostic and restorative instruments, our goal is to remove daily stress for clinicians and improve patient safety. Not only convenient, we have created a product range that is focused on quality and sustainability. Our single use instruments arrive sterile and ready for immediate use. The ergonomic handle is made from ABS plastic and features ridges and grooves for an improved grip. The stainless-steel stem provides additional weight and durability to ensure it is fit for purpose. Precision manufactured stainless steel tips replicate the design and function of re-usable instruments.

Packaging and materials

We are committed to reducing the environmental impact of our products. As such our products and packaging are manufactured with recyclable materials. Additionally, we source recycled cardboard to create our benchtop dispensers and shipping containers.

Quality assurance

All our instrument manufacturing has CE and ISO certification to ensure quality compliance with industry standards.

Tyvek® medical dialysing strip

A Tyvek® medical dialysing strip is incorporated into all our product sleeves. Comprising of flashspun and bonded high-density polyethylene (HDPE) filaments, the strip provides a microbial barrier and exceptional tear and puncture resistance to maintain product integrity during transit.

SINGLE USE

Single Use Dental Instruments Pty Ltd

07 5445 6486 | singleuse.au
sales@sudinstruments.com
1/24 Hoopers Road
Kunda Park, QLD Australia 4556
